

The information in this document contains general descriptions of the technical options available, which do not always have to be present in individual cases and are subject to change without prior notice.

The required features should therefore be specified in each individual case at the time of conclusion of the respective contract.

Local Contact Information

Siemens contact
House number and street name
City and postal code
Country
Phone: +XX-XXX-XXX-XXXX
www.XXXXXXXXXXXXXXXXXX

Global Business Unit

Siemens Audiologische Technik GmbH
Gebbertstrasse 125
91058 Erlangen
Germany
Phone: +49 9131 308 0

Global Siemens Headquarters

Siemens AG
Wittelsbacherplatz 2
80333 Muenchen
Germany

Global Siemens Healthcare Headquarters

Siemens AG
Healthcare Sector
Henkestrasse 127
91052 Erlangen
Germany
Phone: +49 9131 84 0
www.siemens.com/healthcare

Legal Manufacturer

Siemens Audiologische Technik GmbH
Gebbertstrasse 125
DE-91058 Erlangen
Germany

www.siemens.com/hearing



Why miss a decibel
of “Looks like me”?

Siemens Lotus.
For everything you love to hear.

The family of hearing instruments that offers
you great performance and flexibility every day.
www.siemens.com/lotus

Answers for life.

SIEMENS



Welcome to the Siemens Lotus family.

A kind word, a confiding whisper, a mischievous giggle. Decibels of sound full of meaning and emotion. Decibels of sound you don't have to miss – thanks to the powerful, robust and affordable range of hearing instruments from Siemens Lotus.


For hassle-free hearing enjoyment.

Siemens Lotus instruments offer you an array of solutions for mild to profound hearing losses. They are available in a wide variety of Behind-The-Ear (BTE) and In-The-Ear (ITE or "custom") models. And, all Lotus instruments feature advanced feedback cancellation.

Whatever your choice, one thing is clear: Siemens Lotus devices combine affordability with reliability, ease of use and superior sound quality.

New discreet BTE models.

Lotus M is the inconspicuous solution for people with mild to moderate hearing loss. The small Lotus Pro M model can be fitted with an especially thin tube (LifeTube) for greater comfort and discretion. Its stop-end volume control is easy to use. It is also available with directional microphones that make it easier to hear in difficult listening situations.

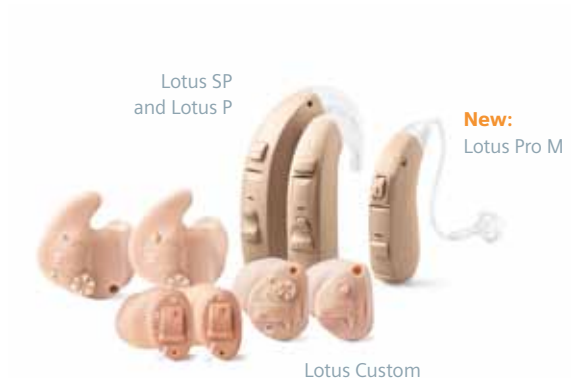

For more confidence in every situation.

Curing a teddy bear's leg is easy – with a bit of care and patience. In comparison, having your hearing instruments repaired can be a nuisance. We would like to spare you the bother altogether. That's why we've made Siemens Lotus products trouble-free and durable. So you'll always feel comfortable and confident when you wear them.

Strong as a bear.

With their robust design, proven components and advanced features, Siemens Lotus instruments have everything you need for highly reliable hearing performance and comfort in the most diverse and demanding situations.

Well protected.

Lotus instruments are resistant to humidity – this protects the interior against damage and minimizes the need for repairs.

Made to last.

Trustworthy performance is the hallmark of the Siemens Lotus products. We pay special attention to details such as the reinforced battery door, the strong metal connector earhooks, volume controls with humidity protection, and reliable push button switches.





Thanks, Doc:
68 decibels



Beautiful:
80 decibels

For easy handling.

Imagination, skill, experience...

We think you should use your talents for the things you really want to do – and not for fiddling around with your hearing instruments. That's why we have put so much work into making them easy to handle. Because real quality and reliability are nothing without ease of use.

Easy from the start.

Siemens provides your Hearing Care Professional with the information and technology necessary to fit your Lotus instruments quickly and accurately.

Convenient details.

Great attention was also given to comfort when designing Siemens Lotus. Clear beep signals are used to indicate low battery as well as volume control and push button adjustments. So you have constant control over the status of your instruments.

Attractive extras.

Some Lotus models are available with a Telephone-coil (T-coil) function to help make phone calls more comfortable. The T-coil prevents the undesirable whistling that can sometimes occur when using the telephone while wearing a hearing instrument.



Easy handling in just three steps.

For discovering the power of small things.

Finding the right instrument for your hearing needs is easy with the Siemens Lotus family. From the smallest Complete-In-the-Canal (CIC) hearing aids to the Super Power BTE models, Lotus delivers the performance you require in a way that is efficient, reliable and discreet.

High output.

The Lotus SP has a high output and gain of 138/80 decibels that is suitable for severe to profound hearing loss. The custom line offers almost invisible In-The-Ear models with more power than ever before. While the attractive directional BTE instruments include Power and Super Power solutions.

Attention to detail.

In order to provide the best sound quality Siemens Lotus instruments use techniques such as feedback cancellation for excellent stability and microphone noise reduction to stop you from hearing the microphones working in quiet environments. The directional microphones also help in difficult listening situations.

Adapting to the environment.

In addition to the “universal” or basic setting, you can select up to three different hearing programs at the touch of a button. There is a program for noisy environments, one for listening to music or TV, and a Telephone-coil (T-coil) program for listening to the telephone.



Lotus P (133/70 dB).
For moderate to severe hearing loss.



Lotus SP (138/80 dB).
For severe to profound hearing loss.



I can't wait:
76 decibels



Family reunion:
64 decibels

For longer hearing enjoyment.

Siemens Lotus instruments are powerful, but energy efficient. They provide great performance at low power consumption to ensure long battery life. By using Siemens batteries, you'll be certain that your instruments remain reliable, safe, versatile and extremely powerful.

Long-life quality.

Even the most powerful and demanding instruments in the Lotus family deliver long-lasting and stable amplification thanks to their economic battery consumption. The excellent battery life is complemented by alert tones when the battery capacity is low, so you won't be caught out without power.



Lotus SP, Lotus P
and Lotus Pro M:
High power performance,
low battery consumption.

For everything you love to hear.

With its robust design, proven components and advanced features, Siemens Lotus has all the elements necessary for highly reliable hearing performance and comfort in the most diverse and demanding situations. Lotus instruments come in programmable and trimmer-operated versions. The Lotus Pro can be programmed using a computer to match your hearing loss characteristics. The Lotus trimmer instruments can be adjusted with controls called trimmers, which are on the instrument itself.

The Lotus family range:

Lotus Pro (Programmable) instruments:

BTE	Lotus Pro M
	Lotus Pro P
	Lotus Pro SP
Custom	Lotus Pro CIC (Complete-In-the-Canal)
	Lotus Pro ITC (In-The-Canal)
	Lotus Pro ITE (In-The-Ear)

Lotus trimmer instruments:


BTE	Lotus M
	Lotus P
	Lotus SP
Custom	Lotus CIC (Complete-In-the-Canal)
	Lotus ITC (In-The-Canal)
	Lotus ITE (In-The-Ear)

The various Lotus instruments differ in power, features and size.


Colors:



Beige

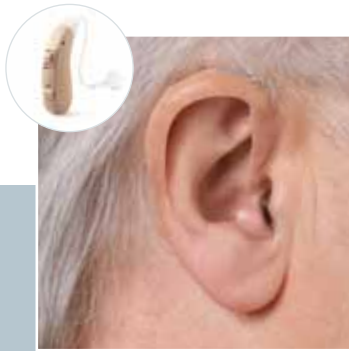

Lotus Pro M with Lifetube*

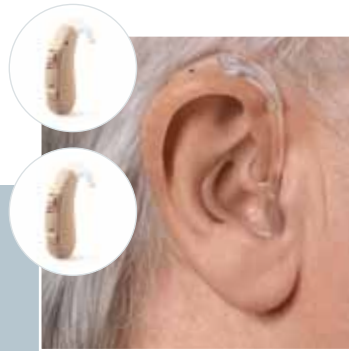

Lotus Pro M with Earhook*
Lotus M**



Lotus Pro P with Earhook*
Lotus P**



Lotus Pro SP*
Lotus SP**



Lotus Pro CIC*
Lotus CIC**



Lotus Pro ITC*
Lotus ITC**



Lotus Pro ITE*
Lotus ITE**



*Lotus programmable instruments **Lotus trimmer instruments

Features and descriptions.



The table below highlights the array of features available in the various Lotus instruments.

	Programmable	Channels	Programs	Directional Microphone	LifeTube	Volume Control	Push Button	T-coil	Audio Input	Antiphase Feedback Cancellation
BTE	Lotus Pro M	2	3+A	✓	optional	Stop-end VC	✓	✓	✓	✓
	Lotus Pro P	2	3+A	–	–	Rocker switch	✓	✓	✓	✓
	Lotus Pro SP	2	3+A	–	–	Rocker switch	✓	✓	✓	✓
Custom	Lotus Pro CIC	2	3	–	–	–	optional	–	–	✓
	Lotus Pro ITC	2	3	–	–	optional	optional	optional	–	✓
	Lotus Pro ITE	2	3	–	–	optional	optional	optional	–	✓
Trimmer										
BTE	Lotus M	2	1+T	–	–	Stop-end VC	✓	optional	✓	✓
	Lotus P	1 or 2	1+T	–	–	Rocker switch	✓	optional	✓	✓
	Lotus SP	1 or 2	1+T	–	–	Rocker switch	✓	optional	✓	✓
Custom	Lotus CIC	1 or 2	1	–	–	Screw Set VC	–	–	–	✓
	Lotus ITC	1 or 2	1	–	–	Stop-end VC	–	optional	–	✓
	Lotus ITE	1 or 2	1	–	–	Stop-end VC	–	optional	–	✓

An elderly couple is shown in a warm, intimate setting. The woman, on the left, has short blonde hair and is wearing a light-colored, textured cardigan over a red top. She is smiling broadly and looking towards the man. The man, on the right, has grey hair and is wearing a blue button-down shirt under a grey cardigan. He is sitting in a wooden chair, holding an open red book, and looking up and to the right with a gentle smile. The background features a framed painting of a landscape with trees and a window with white curtains.

Story time:
54 decibels

Siemens Lotus.
For everything you love to hear.